



eFOLDi User Manual

Model: Power Chair HBLD3-D



SunTech UK Ltd

Serial number:

Please see image showing the location of the serial number.



Contents

1	Welcome	2
2	Terms, conditions and disclaimers	3
3	Product overview	6
4	Components	7
5	Safety guidance	9
6	Folding instructions	14
7	Drive your power chair	18
8	Battery and Charging	23
9	Checks and maintenance	28
10	Refurbishment	29
11	Post-Market Surveillance	30
12	Technical specification	31
13	Appendix	32

Welcome

Thank you for choosing an eFOLDi product as your mobility solution.

I'm Sumi Wang, the founder of SunTech UK and would like to welcome you to the eFOLDi family of users.

SunTech UK manufactures a range of assistive technologies to help you get the most from life. Our mission is to give our customers freedom, friendship and fun with a range of easy to use, convenient and market leading mobility solutions.

Your eFOLDi has been manufactured to the highest standards and quality checked by our own pre-sales technicians before delivery, to provide you with a product that will give many years of freedom and fun with reliable service.

Please take some time to read the user manual and the supplied documents before using your eFOLDi for the first time, this will help you understand the product and help you obtain maximum enjoyment and safety from your new eFOLDi product(s). You will quickly become familiar with the features of your eFOLDi with a little practice.

SunTech UK prides itself on its high-quality products and customer support so please remember to register your product(s) at: <https://efoldi.com/register-warranty>

This will help us to assist you in the event of any difficulties you may experience in the future. I am sure you will find your eFOLDi to be a reliable friend and companion and I thank you for your custom and your trust in our products.

Yours,



Sumi Wang

Terms, conditions and disclaimers

The terms 'us', 'we', 'our', 'SunTech UK' and 'SunTech UK Ltd' mean SunTech UK Limited, which is a UK registered company, limited by guarantee, registered with Companies House, number 06906908. The term 'website' means the website found online at www.efoldi.com. The terms 'you', 'your', 'user' and 'rider', mean the person using the eFOLDi notwithstanding any person(s) handling the eFOLDi on behalf of the user, if different to the user, including general handling, maintenance, repairs and so on.

Purchase agreement

By accepting delivery of this eFOLDi, you undertake to not change, alter or modify this eFOLDi and/or any parts supplied with it or subsequently for it, and/or render inoperable or unsafe, any guards, shields or other safety features of this product; fail, refuse or neglect to install any retrofit kits from time to time, as may be provided by SunTech UK Ltd to enhance and/or preserve the safe use of this product. You also agree to accept and abide by/with the terms and conditions contained herein and as may be revised from time to time and updated on our website, being immediately applicable. By using the eFOLDi and/or accessories and/or components supplied with/or for the eFOLDi, the user acknowledges having read and understood the use(s) and limitations of the product(s) without exceptions. This shall not act as a waiver of any implied warranties, which the customer retains, as provided by law.

Warranties will be void and SunTech UK Ltd will accept no responsibility whatsoever or howsoever, in instances of misuse and/or abuse of the products, in part or in full, by the user and/or any unauthorised persons, which may include (but not to be limited to); the user(s) failing to comply with all legal requirements and guidelines, wherever and/or however it is used or intended to be used, lack of/or inappropriate maintenance, maintenance (apart from the routine maintenance as described in this manual by the user) by anyone not approved to do so, unapproved repairs, unauthorised disassembly of the vehicle and/or components includ-

ing (but not limited to) battery pack(s) and controllers, unapproved modifications or alterations, failure to ensure correct folding and unfolding procedures are followed precisely and in full, damage or injury howsoever caused, failure to apply all reasonable and appropriate safety precautions at all times, improper use of the vehicle and/or parts and/or accessories being used by more than one person at any time, use in any competitive sport, racing, stunts, jumping or other similar activity.

Our liability and the liability of our suppliers, to you/or any third parties in any circumstances is limited to the user's cost of qualifying products purchased. Qualifying products as acceptable for the return, as outlined in the return policy, as stated on the website online at: www.efoldi.com/return-policy.

Information in this user manual is subject to change without prior notice. For the latest updates and information, please go to our website online at: www.efoldi.com.



The power chair should NOT be used by people:

- A person who has an intellectual or developmental disabilities;
- With slow reaction and who have difficulty in handling power mobility devices;
- Who cannot operate all controls fully and as intended; *
- Who cannot steer competently and appropriately at all times and in all conditions;
- Who cannot maintain appropriate awareness of all relevant circumstances for each journey and take appropriate action as may be required to avoid mishap;
- Who cannot lift it in whole or in part as may be required dependent upon the user's needs*.

* Unless accompanied by a person sufficiently able and competent to carry out that task.

Product compliance:

Safety:	EN12184
EMC:	ISO7176-21
Risk management:	ISO14971
The European Medical Devices:	Regulation (EU) 2017/745
According to EN12184 Class:	Class A

Shipping and delivery:

By checking against your sales order or invoice, make sure your delivery is complete, as some components may be individually packaged. If you do not receive a complete delivery, please contact us immediately. Where damage has occurred during transport, either to the packaging or the contents, please contact the company responsible.

Product overview

Thank you for purchasing our eFOLDi power chair. This ultra-light electric wheelchair has been carefully researched and developed to be light, energy efficient and easy to operate. Please read this manual carefully before use, so that you are familiar with all the various functions, and understand the required maintenance. This will ensure that the power chair remains in good condition and is safe to use.



(a) Power chair mode




(b) Storage mode

Model: Power chair HBLD3-D

Intended for use:

This Class A Power Chair is intended for use by disabled persons and/or individuals with mobility impairments to address walking restrictions. It is designed for indoor and limited outdoor use on smooth, level surfaces - such as:
Homes, apartments, and assisted living facilities;
Shopping malls, offices, and hospitals.

Components

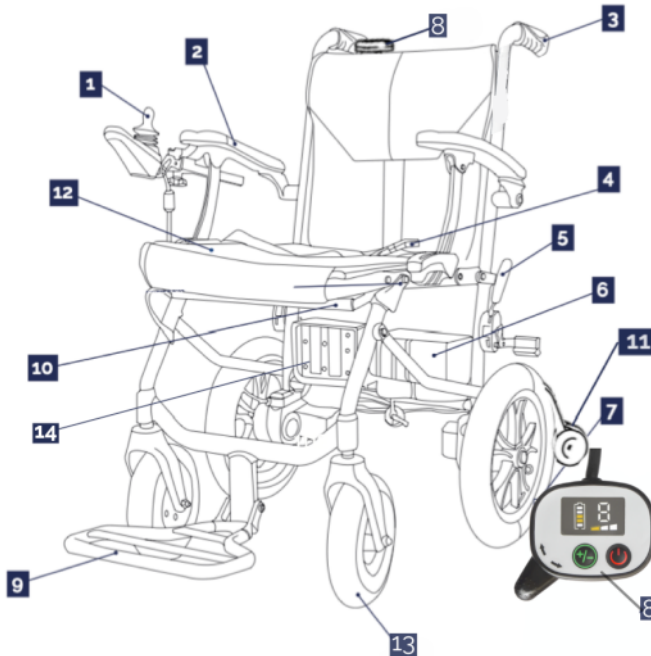

1	Joystick control	2	Armrest
3	Push handle	4	Seatbelt
5	Handbrakes (both sides)	6	Batteries
7	Rear wheel 12 inch	8	Attendant control
9	Footrest	10	Folding mechanism bar
11	Anti-tip wheels	12	Seat cover
13	Front wheel	14	Control box

Table 1: Model HBLD3-D Components

Control components



1	Joystick control	2	Speed adjustment (slower)
3	Battery power display	4	Power switch
5	Horn	6	Speed adjustment (faster)
7	Speed indicator	8*	Attendant control
9	Reversing sound on/off		

Table 2: Power chair HBLD3-D controls

* Attendant control only works when the main power switch (4) is turned on. Push any buttons on the attendant control to sync it with the Joystick before moving off.

Max speed**	Joystick	Attendant control
1	2 km/h	2 km/h
2	3 km/h	3 km/h
3	4 km/h	4 km/h
4	5 km/h	NA
5	6 km/h	NA

Table 3: Power chair speed settings

** Speed may vary due to manufacturing tolerances and usage conditions.

Safety guidance



Attention

- It is the user's responsibility to read the manual and supplied documents before using the eFOLDi and to follow all of the recommendations contained within, to fully comply with all local laws of their countries and localities and to use the eFOLDi with due diligence and regard to other road users and pedestrians.
- The eFOLDi is simple to use, but its folding mechanism can cause significant injury if used inappropriately and handled without due care and diligence. To avoid any risk of injury ensure hands, fingers, jewellery, clothing and other objects are not allowed to interfere with the folding mechanism or moving parts of the eFOLDi and study the manual carefully before you start operating your eFOLDi.
- It is advised that the user familiarise themselves with the folding and unfolding procedure before venturing outside. The user should also ensure that they familiarise themselves with the controls and handling of the power chair in a safe area away from obstacles and other pedestrians or road users.
- Fully charge the battery(ies) before first use and follow the guidelines in this manual for prolonging the life of your battery.



Do not get on or off the eFOLDi power chair when the power is on, or when the power chair is in freewheel mode or with the handbrakes disengaged, or there is no carer to provide assistance.

Please do NOT riding over rough, cobbled, muddy or loose gravel surfaces. Your eFOLDi is designed for smooth, flat surfaces and may struggle with uneven pavement, gravel, or grass.

The stopping distance on slopes can be significantly greater than on level ground; The power chair surface temperatures can increase when exposed to external sources of heat (e.g. sunlight)

Cannot be driven on the road.

Be sure to slow down before turning. Do not make sharp turns at high speed. Do not park the eFOLDi power chair on a ramp. Never make turns on a slope.

Any vehicle can be dangerous to its rider/driver and other road users and pedestrians. Always ride/drive responsibly, respectfully and with due care to others.

Your eFOLDi power chair is NOT a toy and is not recommended for people under the age of 14. Do not allow children to play with or on the power chair. The internal folding mechanism can cause injury to small hands.

Do not operate your eFOLDi power chair until you have read and fully understood this manual.

Do not operate your eFOLDi power chair until assembly and inspection is complete. A person who has an intellectual or developmental disabilities, slow reactions and difficulty in handling powered devices should not use the eFOLDi power chair on their own.

Do not drive your eFOLDi power chair in violation of national and local traffic regulations.

Do NOT rush when folding or unfolding the eFOLDi power chair.

Do not use the eFOLDi power chair on sandy or soft ground. Avoid driving on slopes greater than 6 degrees or over obstacles greater than 4 cm in height.

Do NOT disassemble or change parts of the eFOLDi power chair or replace parts not manufactured by the manufacturer without the manufacturer's permission.

Never operate the power chair without the anti-tip wheels securely in place.

Do NOT for any reason put your hands into the mainframe of the power chair. This is because your hands and fingers can get trapped in the levers and pivots and serious injuries can be caused.

The use of accessories, transducers and cables other than those specified or provided by the manufacturer of this equipment could result in increased electromagnetic emissions or decreased electromagnetic immunity of this equipment and result in improper operation.

Do not use the power chair on a ramp/slope if the battery is out of power without applying manual brakes, it does not have any electric brake without power.

If there is NO charge within the batteries (or the batteries have been removed from the Power Chair), the Power Chair will free wheel at any time, provided the manual handbrakes are released.

The Attendant Control has no impact on the braking system of the Power Chair.



When unfolding/folding your power chair ensure:

- The eFOLDi is positioned on a level floor or ground.
- Main switch is switched off.
- All latch pins are properly secured in place, as stated on the instructions and shown on the illustrations.
- Handbrake is disengaged.



Safety check before use:

It is highly recommended that you perform a safety check before each use to make sure your eFOLDi operates smoothly and safely.

- This power chair is not designed for use as a seat in motor vehicles. Always transfer to a proper vehicle seat and secure the chair separately.
- Check for correct tyre inflation. Maintain but do not exceed the manufacturer recommended tyre pressure rating indicated on the pneumatic tyres;
- Check all electrical connections. Make sure they are secure and not corroded;
- Check and test the brake system;
- Check frequently for any loose fasteners, fixings and other components.
- Check Anti-tip wheels are securely attached.



Storage conditions:

Your Powered Wheelchair should be stored in a dry place, free from temperature extremes. When storing, disconnect the battery from the Powered Wheelchair. See "Batteries and Charging."

WARNING! If you fail to store the unit properly, the frame can corrode and the electronics can be damaged.

Batteries that are regularly and deeply discharged, infrequently charged, stored in extreme temperatures, or stored without a full charge may be permanently damaged, causing unreliable performance and limited service life. It is recommended that you charge the battery periodically, every 3 months, throughout periods of prolonged storage to ensure proper performance.

- Storage temperature: -25°C -to 50°C
Storage relative humidity: 60 % ± 20 %
- Some of the parts of the power chair are susceptible to change in temperature. The controller can only operate in temperature that ranges between - 25°C to 50°C.



EMC warning statement

The interference from radio wave sources such as radio stations, TV stations, two way radios or cellular phones can cause the mobility product to release its brakes, move by itself, or move in unintended directions. It can also permanently damage the mobility power chair control system.

The intensity of the interfering EM energy can be measured in volts per meter (V/m). Each mobility product can resist EMI up to a certain intensity. This is called the "EMI immunity level". The higher the EMI immunity level, the greater the protection from EMI. Currently the eFOLDi power chair has an immunity level of at least 20V/m, which provides useful protection from the common sources of EMI.

The sources of radiated EMI can be broadly classified into three types: 1) Hand-held portable transceivers. (normally with antenna mounted on the device) 2) Medium-range mobile transceivers, such as those used in police cars, fire trucks, ambulances and taxis. (normally with antenna mounted on the vehicle) 3) Long-range transmitters and transceivers, such as commercial broadcast transmitters (radio and TV broadcast antenna towers) and amateur (HAM) radios.

Other types of hand-held devices, such as cordless phones, laptop computers, AM/FM radios, TV sets and small appliances such as electric shavers and hair dryers are usually not likely to cause EMI problems.

Use of this electric wheelchair adjacent to or stacked with other electrical equipment should be avoided because it could result in improper operation. If such use is necessary, this equipment and the other equipment should be observed to verify that they are operating normally.

Portable RF communications equipment (including peripherals such as antenna cables and external antennas) should be used no closer than

30 cm (12 inches) to any part of the electric wheelchair, including cables specified by the manufacturer. Otherwise, degradation of the performance of this equipment could result.

If the use location is near (e.g. less than 1.5 km from) AM, FM or TV broadcast antennas, before using this equipment, it should be observed to verify that it is operating normally to assure that the equipment remains safe with regard to electromagnetic disturbances throughout the expected service life.



MR Safety statement




"The product is MR Unsafe. The device presents a projectile hazard."

Folding instructions

To a Power chair mode




Please study the instructions and image illumination to familiarise yourself with the folding procedures and please only operate one moving part at the same time.

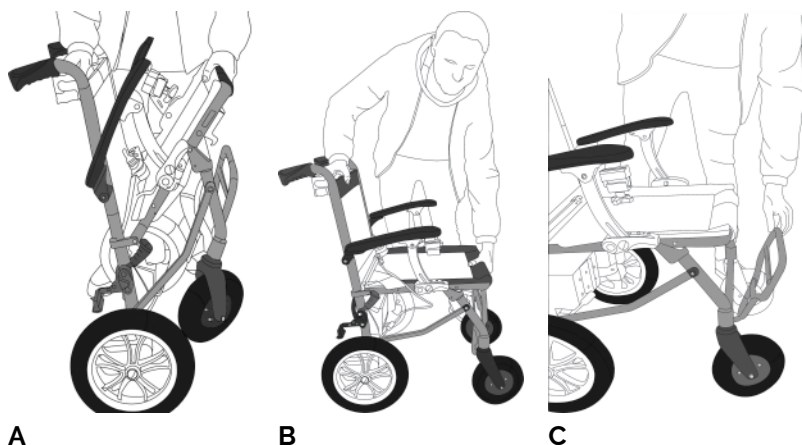


A Once removed from the packaging place the wheelchair on a stable surface. Place one hand on the top of the backrest and the other hand on the folding mechanism bar underneath the seat, in the front part of the folded power chair.

B Finally push your two hands apart until the seat is open at 90° degrees and you hear a distinct "click" sound.

C Unfold the footrest, from a vertical position to a 90° degrees horizontal position.

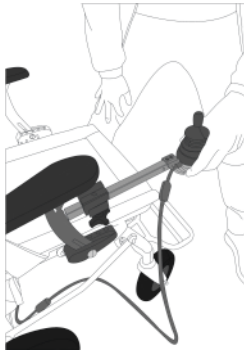

D

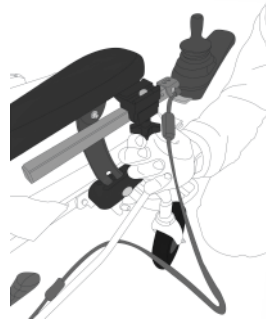

E

D Slot the joystick arm inside the frame and push it into a comfortable horizontal position, where it is easy for your hand to reach the joystick.

E Once you are happy with the horizontal position of the joystick arm, lock it by tightening the clamp. Finally attach Anti-tip wheels securely on the back of the power chair before driving off.



Power chair mode


To a Storage mode




Please study the instructions and image illumination to familiarise yourself with the folding procedures and please only operate one moving part at the same time.



A

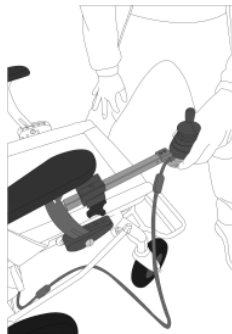

B

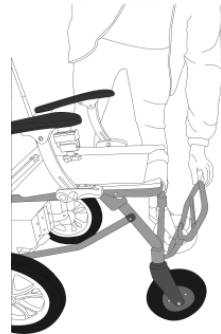

C


A Loosen the tightening clamp, to release the joystick arm.

B Take out the joystick arm and place it on the power chair seat.

C Fold the footrest, at 90° degrees position vertically.



E

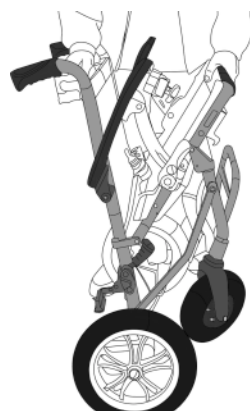

F

E Then place one hand at the top on the backrest bar and the other on the folding mechanism bar underneath the seat.

F By pressing the folding mechanism bar and bringing your hand together to close and fold the eFOLDi Power chair. Your power chair is ready to be stored or transported.



Front wheel should face backwards for proper folding.



Storage mode

Drive your power chair



Safety check before use - see Safety Guidance

eFOLDi power chair is suitable to use indoors and is designed to drive on smooth even surfaces. When driving your eFOLDi, choose a route with good road conditions. Unevenly paved roads or paths, potholes and other hazards will reduce the life of the power chair and impact on your safety.

Any deviation from level/horizontal surfaces should not exceed a gradient of 1:20 (3 degrees). For any climbing there are many factors that need to be considered. Such as: the tyre pressure and condition, the surface condition (if wet, icy or slippery, and this includes leaf cover and/or other debris, we advise against driving the slope), battery charge and the weight of the user, all of these will affect performance.



Caution:

The stopping distance on slopes can be significantly greater than on level ground; The power chair surface temperatures can increase when exposed to external sources of heat (e.g. sunlight)

ALWAYS APPLY MANUAL BRAKES BEFORE GETTING IN TO OR OUT OF THE POWER CHAIR

Joystick controller



**Perform the safety check before each use
- see Safety Guidance section.**

Operate the power chair

- To Start - Switch on the main power switch (4) and press the power button on the joystick to switch on the power chair.
- To Stop - Switch off the power switch (4) on the joystick and apply manual brakes, which are located on both sides of the chair. When parking on a ramp/slope, the electric brake may not always be strong enough. Always use the manual brake in conjunction with the electric brake.
- Adjust the speed - Use speed adjustment button to reduce (2) or increase (6) the speed.
- Drive forward - Push the Joystick controller stick (1) forward.

- Turn left/right - Push the Joystick controller stick (**1**) left/right.
- Drive backward - Push the Joystick controller stick (**1**) backward.
- Use attendant control - Power on the attendant control (**8**) and push the lever forward to go forward, push back the lever to go backward.



When Attendant Control (**8**) is in use, the chair no longer responds to the Joystick controller (**4**).

- Joystick controller (**1**) controls the direction and speed at the same time, slowly push the joystick towards the driving direction, then the eFOLDi Power chair will start to move, push the joystick further and you will notice greater acceleration.
- Speed adjustment - slower (**2**) It is positioned just above the control handle on the left hand side, and it decreases the speed (from 5 to 1) to make the eFOLDi power chair easier to control.
- Power display (**3**) Shows the battery charge remaining.
- Power button (**4**) Switches the power on and off by pressing the button (Make sure the main power switch is on).
- Horn (**5**) Positioned right in the middle of the control panel, use this to warn other pedestrians of your presence.
- Speed adjustment - faster (**6**) It is positioned just above the control handle on the right hand side, and by pressing the button it increases the speed level (from 1 to 5) of the power chair to go faster. Always be aware of your speed.
- Speed display (**7**) Indicates what speed is currently selected, between 1 to 5.

- Attendant control (8) Allows an assistant to power the chair for easy handling. Also it has its own power on/off button.
- Reversing sound on/off button (9) Allows you to turn the reversing beep on or off.

Freewheel/motor bypass mode

Switch off the power switch from joystick control on the flat ground for 2 seconds or lift the rear wheels off the ground for a couple of seconds, it will give you freewheel/motor bypass mode, but you will still feel some resistance.

Braking mode



Do not use the power chair on a ramp/slope if the battery is out of power without applying manual brakes, it does not have any electric brake without power.

If there is NO charge within the batteries (or the batteries have been removed from the Power Chair), the Power Chair will free wheel at any time, provided the manual handbrakes are released.

The Attendant Control has no impact on the braking system of the Power Chair.



Manual handbrakes

When on a slope, when coming to a stop. Do NOT turn the power

off until the manual brakes are applied. If the power is turned off prior to applying the manual brakes, the electric brake will hold the chair up for to 15 minutes (provided there is charge within the batteries). After 15 minutes, the Power Chair will start to roll forward, though there will be some friction slowing this movement provided there is still some charge within the batteries. Apply manual brakes immediately to stop the forward movement.



ALWAYS APPLY MANUAL BRAKES WHEN THE POWER CHAIR IS STATIONARY, ESPECIALLY IF ON A SLOPE, AND ALWAYS BEFORE GETTING IN TO OR OUT OF THE POWER CHAIR

Transporting a power chair:

To safely transport the power chair, first ensure it is powered off. Fold the chair according to the manufacturer's instructions, securing any loose parts. Lift using the designated handles or frame sections, avoiding delicate components. Place it in the vehicle with the battery properly secured. If necessary, use straps to prevent movement during transport.



Caution:

The product is not intended to be used as a seat in a motor vehicle. You will need to take a seat on your mode of transport.



Battery and Charging



Please read the instruction before using the charger. The eFOLDi power chair battery has been designed and built specifically for the eFOLDi power chair. This battery must not be used for any purpose other than as stated in this manual. Any misuse or dismantling of this battery may damage it physically and/or electrically, internally and/or externally and will not be covered by warranty. The battery contains no user serviceable or replaceable components. The eFOLDi power chair battery is actually a battery pack, comprising a number of lithium ion cells, but merely termed 'battery' for ease throughout this manual. This battery is different in characteristics and charging requirements to the old and still quite common lead/acid or lead/gel type. Never use a lead type charger to charge it.

It is highly recommended to ensure the charger is placed only on heat resistant surfaces and nothing which could contribute to a fire (e.g. carpet, floorboards, wood, laminate or plastic worktops etc)- the charger will heat up during a normal charging cycle, as will the battery.

To ensure best long-term use of the battery, please observe the following safety rules at all times.

- Connection dock must be handled with care to avoid damage.
- Charge the battery ONLY with the approved charger.
- Never use a charger for lead type batteries on the eFOLDi battery.
- Never place the battery near, on or against any heat source.
- Ensure no metal parts can touch (or 'short') battery or charger terminals, connectors, pins or any other electrical conductor(s) or conducting parts.


- Avoid excessive physical impact and vibration of the battery.
- Do not attempt to dismantle, dissect or deform the battery(ies) or case(s).
- Do not immerse the battery in any liquid (including water).
- Do not mix different types of batteries.
- Keep the battery out of reach of children.
- Do not use modified or damaged chargers.
- Do not leave the battery on charge overnight.
- Store the battery in a cool, dry, well ventilated area.
- For disposal of the battery must be in accordance with local laws and regulations, please contact your local dealer or SunTech UK Ltd (legal manufacturer) for advice.
- Batteries should not be stored for long periods without charging. It is recommended the battery to be charged monthly to prolong its life.
- Never store the battery when it is significantly discharged.
- For the benefit of international travellers, your approved eFOLDi charger can be used with any mains voltage supply from 100 to 240 Volts, with a suitable travel adaptor.

Battery power display

Switch on the power, the power display light in the front row will be on, full power will show 5 lights (2 red, 1 yellow, 2 green). When the two red lights are flashing, this indicates that the battery requires charging.



To remove the battery for charging, push the red button down, hold the battery on the outside, and release from the Power chair. To put it back, slide the battery into the battery guide slot, push it until you hear a click, then the battery is inserted securely.

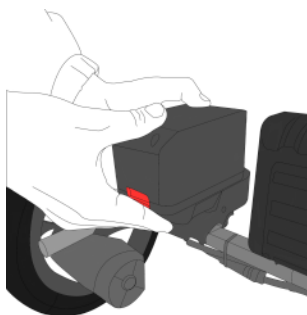

Caution:

Please charge the battery once a month if not in use (Power chair batteries ideally should not be stored in the power chair).

Do not dispose of batteries in regular trash — recycle via designated facilities or hazardous waste programs.

Battery dimensions: 12cmx17.5cmx9cm


Battery discharge capacity: (501st cycles) $\geq 70\%$

Charger and Charging

Use only the charger originally supplied with your new eFOLDi power chair or replacement charger supplied by SunTech UK Ltd for the use with your product. The use of a non-approved charger will invalidate any warranty and may cause serious damage to the battery and/or your eFOLDi product. The charger should not be used with/or for other devices as damage caused to the charger and/or any other device by doing so will not be covered by our warranty.

Do not leave the battery on charge for more than 24 hours, as this will reduce the life and performance of the battery.

See the image below. Insert the charger into the charging port to start charging.

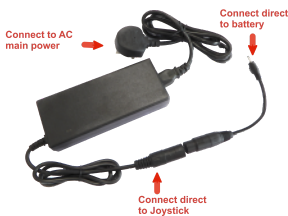

You can charge the battery in 2 ways:

1. Charge via Joystick - plug the charger lead into the base of the joystick controller and connect the AC plug to a power outlet. This will charge both batteries simultaneously if both batteries are on the bays.
2. Charge directly to battery - battery may be charged separately after removing it from the chair by using the short 'pigtail' adaptor which needs to be connected to the charger lead first.



Note: see the back of the charger for the charging status.



Stop Using a damaged battery or battery charger immediately if your charger has the following signs:

1. Frayed/Exposed wires
2. Burn marks, melted plastic, or unusual odors
3. Intermittent charging or overheating
4. Cracked casing or loose connectors

Continuing to use a faulty or damaged charger can lead to severe risks, including:

1. Fire or Explosion:
2. Electric Shock:
3. Battery damage.
4. Device damage.

Checks and maintenance



- Perform the following inspections prior to using your eFOLDi daily: For pneumatic tyres, make sure you check the correct tyre inflation and maintain without exceeding the P.S.I./BAR/kPa air pressure rating indicated on each pneumatic tyre (note: the front and rear tyres may have different pressures). Also be sure to use the correct pressure scale and according number on your inflator/pressure tester.
- Perform the following inspections weekly: Check the fastening condition of all nuts, bolts and screws to make sure all the parts are properly secured and are not missing or loose. Check all electrical connections. Make sure they are tight and not corroded. Electromagnetic brake - is used as a driving brake. The braking distance should be less than 1.5 metres after you release the joystick to stop. Tyre pressure - See your tyres for this information.
- Cleaning and maintenance Periodic cleaning of your eFOLDi will help to maintain its appearance and will both extend the life of the vehicle and improve resale value. Metal and plastic parts can be wiped with a damp cloth using a mild detergent.
- Have your eFOLDi regularly serviced each year.

Expected Service Life: The Power Chair has an expected service life of 5 to 7 years under normal usage conditions. The actual lifespan may vary depending on factors such as frequency of use, maintenance, and environmental conditions.



Caution:

Oil and other lubricants must NOT be used on the following parts: Brake pads, brake cable, tyres.

Refurbishment



- Refurbishment for re-use

The product is suitable for refurbishment and re-use.

Actions to be carried out:

Cleaning and disinfection. Please see section "Checks and maintenance".

Inspection according to service plan.

SunTech UK is the sole entity responsible for the refurbishment and reuse of devices under its brand, in compliance with applicable regulations and sustainability commitments.

Post-Market Surveillance

As part of our commitment to product safety and compliance, we conduct post-market surveillance (PMS) to monitor the performance and reliability of this power chair after sale.

Reporting Issues and Safety Concerns:

If you experience any issues related to performance, durability, battery life, or safety, please contact us immediately:

Customer Support Hotline:

+44 (0) 20 3143 5168

Email:

services@efoldi.com


Website:

www.efoldi.com

Address:

**25 Ormside Way,
Redhill, Surrey, RH1 2LW**

Technical specification

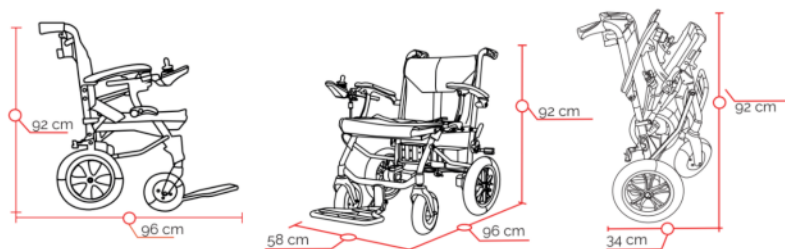

Frame material	Aluminium Alloy
Dimension	
Packaging size (L/W/H):	37x60x94cm (15x24x37")
Folded size (L/W/H):	34x58x92cm (13x23x36")
Unfolded size (L/W/H):	96x58x92cm (38x23x36")
Front wheel size:	20cm (8")
Rear wheel size:	30cm (12")
Seat height:	48cm (19")
Ground clearance:	4cm (1.5")
Weight	
Weight with battery:	17kg(38lbs)
Weight without battery:	15kg(33lbs)
Performance	
Max speed:	≤6kmh
Maximum range capacity:	9miles (15km)
Maximum carry capacity:	135kg (297lbs/21stone)
Maximum Incline:	1:20 Gradient (3 degrees)
Turning radius left/right:	102 (40")
Maximum kerb height:	3cm (1.2")
Components capability	
Motor:	2 x 200W, 24VDC
Battery:	2 x 5.2Ah 24VDC Lithium Ion
Charger:	Input 110-240VAC, 50/60Hz, output 29.4V 2A
Brake:	Electric brake and manual brakes
Power switch:	Power button on joystick
Number of Wheels:	4
Tyre type:	Solid

Table 4: Specification-Power chair HBLD3-D

Legal Disclaimer Notice:

Specifications may vary due to manufacturing tolerances and variations in use conditions.

Appendix



SunTech UK Ltd


25 Ormside Way, Holmethorpe Industrial Estate, Redhill, Surrey, RH1 2LW

TEL: +44 (0) 2031435168



MedPath GmbH

Mies-van-der-Rohe-Strasse 8, 80807 Munich, Germany

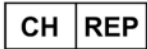


ALBO-Healthcare GmbH

Alte Steinhauserstrasse 19, CH-6330 Cham

Symbols used on the label

	European CE Marking		Manufacturer
	Authorized representative in the European community		Medical device
	Serial number		Unique Device Identification
	Date of manufacture		Refer to instruction manual
	Swiss authorised representative		UK Conformity Assessed

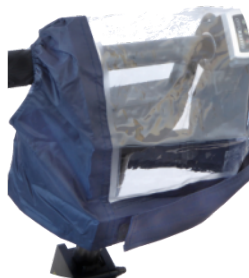

(eFOLDi accessories)
www.efoldi.com/shop



SunTech UK Ltd
 25 Ormside Way, Redhill,
 Surrey, RH1 2LW, UK

**UK
CA**



MedPath GmbH
 Mies-van-der-Rohe-Strasse 8
 80807 Munich, Germany

CE



ALBO-Healthcare GmbH
 Alte Steinhauserstrasse 19
 CH-6330 Cham

Official website address: www.efoldi.com
 Service email: services@efoldi.com
 Facebook: www.facebook.com/efoldi.co.uk