

eFOLDi[®] Navigator Control Box Change

STEP 1




After unplugging and removing the battery, lift the control box out of the elasticated holder and unplug the attendant control cable.

STEP 2



Unplug the joystick cable.

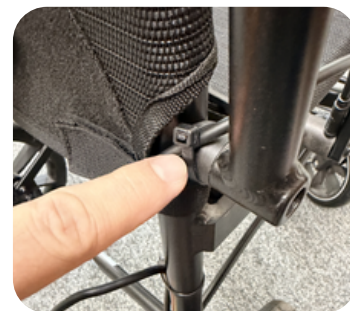

The cable tie securing the cable to the frame will need to be cut.

STEP 3



Locate the two motor to control box cables and unplug them.



Step 4



The faulty control box can now be removed, and the replacement installed.

Step 5 – Ready to Ride



Dealer Portal efoldi.com/dealer-pack

We've made it quicker and easier to access everything you need to sell, demonstrate and support the eFOLDi range.

Our newly redesigned Dealer Portal has been organised around the way you work, with resources grouped by product so you can quickly find exactly what you need - whether you're helping a customer choose the right model, training your team or planning your next marketing campaign.

What's new?

Filter by Product

Find everything you need for each eFOLDi model in just a few clicks.

Enhanced Asset Library

Product images, lifestyle photography, CGI, videos and more—all organised by product.

Marketing Toolkit

Ready-to-use social media posts, guides and promotional materials.

Sales & Support Resources

Brochures, manuals, Quick Start Guides, servicing checklists and tech spec comparisons.

Travel & Accessibility Hub

Air travel forms, checklists, public transport guidance and accessible travel videos.

Consumer Product Profiles

Help customers find the eFOLDi that's right for their lifestyle.

Export Resources

Translated marketing materials and social media templates for international partners.

For assistance, please contact eFOLDi's Customer Care team on +44(0) 203 1435168 or email the team at hello@efoldi.com

