



eFOLDi Cliq

Lightweight mobility solutions

For wherever life takes you

www.efoldi.com



Scan for more
information

**Wherever life takes you,
we'll be right here, making
sure you have the mobility
and freedom to enjoy
every moment.**

**As a recognised leader in the design and
manufacture of portable scooters and powered
wheelchairs, eFOLDi strives to help people
rediscover their freedom and independence.**

eFOLDi's award-winning mobility solutions are not only easy to store and transport, they help make travel simple and achievable – whether that's by air, land or sea.

Innovation

Our products are manufactured using the highest quality materials to ensure they're robust enough to withstand the toughest of conditions and terrain. Our frames are made from high-grade aluminium or magnesium alloys, our high-specification brushless motors deliver exceptional performance, and our batteries use advanced lithium-ion technology, so they're not only durable, they offer outstanding levels of reliability.

COMPASS



CLIQ



LITE



EXPLORER



NAVIGATOR



Why eFOLDi stands apart



Lightweight

Carry and lift
with ease




Warranty

Backed by full
support & warranty



Compact

Portable for greater
independence

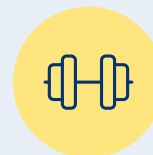

Robust

Peace of mind with
strong, durable
materials



Comfort and portability, clicked together.

The eFOLDi Cliq has four parts each under 10 kg, designed with EZ-Cliq connections for quick, independent assembly.


Enjoy longer trips with a smooth, supportive ride.

For those who value stability, a comfortable seat, and a planted wheelbase, with the flexibility to take the scooter apart for easy transport.



Lightweight

31kg total weight with 4 easy-lift parts

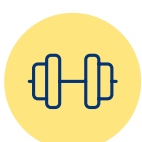

Robust Build

Carries users up to 120kg/265lbs/19 Stone



Range

up to 9miles/14km



Warranty

Backed by full support & warranty*

*terms and conditions apply



Ideal For: individuals with moderate mobility issues who need extra support and stability, seeking a portable scooter that offers increased confidence, security, and ease of use for getting around.



Your Confident Companion on Every Journey

Designed for comfort, stability, and everyday independence.

Extra comfort as standard

Flip-up armrests and practical storage solutions for real-world convenience.

Speed settings

Tortoise/hare mode lets you choose slower precision or faster travel

WigWag controls

Intuitive dual-lever system for simple, confident operation

Delta-shaped handlebars

Ergonomic design for secure grip and comfortable steering

Comfort-focussed seat

Swivel for easy transfers, height adjustable for your ideal position

Robust

Carries users up to 120kg/265lbs/19 Stone

Front and rear lights

For extra visibility

Suspension

Front and rear suspension absorbs bumps and uneven terrain for a smoother, more comfortable ride.

Spacious Foot Area

Helps prevent cramping & discomfort.

Automatic braking

Slows upon throttle release, ideal for users with slow reaction times.

Automotive grade Lithium-ion battery

With 9 mile / 14km range
Removeable for on or offboard charging

Solid tyres

With no pumping requirement & multi-terrain suitability

Handling

Individual parts weigh <10kg with easy-Click connections for effortless assembly and transport



Personalisation

Interchangeable coloured accents as standard to customise your ride



"With the eFOLDi I'm able to visit places I haven't been able to get to for a while. I'm very happy to tell people – buy an eFOLDi."

Peter, Verified Buyer 

Clicks Together

- Modular design: easy assembly and transportation
- Lightweight parts: <10kg with Easy-Cliq connections
- Snap-together simplicity: hassle-free setup
- Fits everywhere: car and home storage without stress

Clicks With You

- Comfort-first seating: adjustable, supportive seat with all-round suspension
- Confident Controls: WigWag levers, ergonomic delta handlebars, simple speed selection
- Stable handling: secure four-wheel chassis for smooth everyday travel
- Made to suit you: personalise with accessories and colour accents

Personalisation

with all interchangeable accent options included as standard, you can mix and match to suit your style and season.



eFoldi Yellow



Moonlight Silver



Harbour Blue



Cliq Specifications*

Maximum Speed	4mph / 6kph
Maximum Range	9miles / 15km
Maximum User Weight	19 Stone / 120kg / 265lbs
Climbing Ability	1:10 Gradient
Turning Radius	121cm (47.6")
Unfolded Dimensions	100 x 61 x 88cm 39 x 24 x 34.6 inches
Individual Part Weight	3.6kg/7.9lbs 9.3kg/20.5lbs 9.8kg/21.6lbs 7.3kg/15.9lbs
Combined Weight (excl. battery)	28kg/61lbs
Combined Weight (incl. battery)	31kg/69lbs
Type of Use	Pavement, Indoor
Braking System	Electro Magnetic
Airline friendly	Yes

*Specifications may vary due to manufacturing tolerances



7 Year Battery + 2 Year Scooter Warranty***

***terms and conditions apply



Scan for more
information

eFOLDi 

For more information...

Call us today on
0800 999 8123

Or email us at:
sales@efoldi.com

Your local eFOLDi retailer:



SunTech UK Limited t/a eFOLDi, Unit H The Courtyard
Business Centre, Lonesome Lane, Reigate, Surrey, RH2 7QT